



EXPERTS IN PATIENT HANDLING SOLUTIONS



CAR DOOR PULL STRAP

This guide covers the product DS0540C.


This **Car Door Pull Strap** is for those who struggle to shut their car doors because they can't reach. Suitable only for car doors with an open grip so ensure that it's compatible with your car door handle.

USER INSTRUCTIONS



1

Loop the **Car Door Pull Strap** through the handle of the car door as shown, holding the **Car Door Pull Strap** tight.

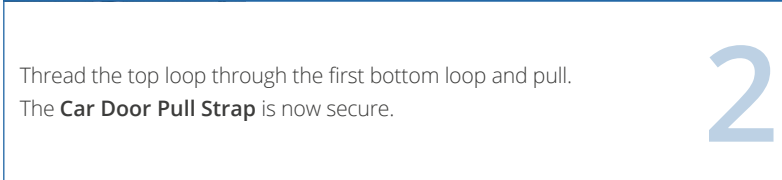


2

Thread the top loop through the first bottom loop and pull. The **Car Door Pull Strap** is now secure.



3

Sitting in the car, you can now use the **Car Door Pull Strap** to pull and close the door.

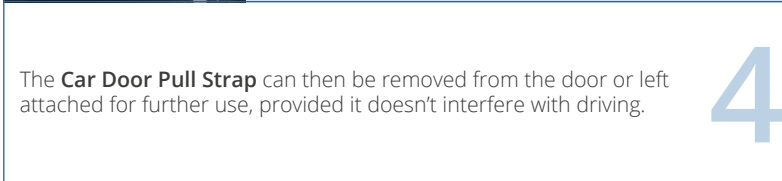

4


The **Car Door Pull Strap** can then be removed from the door or left attached for further use, provided it doesn't interfere with driving.

FABRIC CONTENT

100% Polypropylene


STORAGE AND HANDLING

- Store in a clean dry environment that is not subject to excessive heat.
- Keep away from sunlight.
- Store flat and avoid placing heavy objects on top of the product.

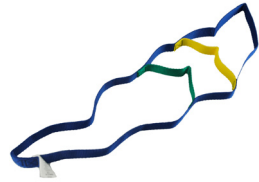

WARNINGS AND PRECAUTIONS

- Do not leave or drop on the floor; it is a slip or trip hazard.
- Test all slide sheets for slippiness before use. If over-laundered and they don't slide please discard and replace.
- Inspect all products, including handles, before use for signs of damage or tears..

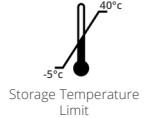

CLEANING INSTRUCTIONS

- Wash as per laundry instructions on the previous pages.
- Can be wiped clean, between use with the same patient, with natural water and detergent or wipes that don't contain alcohol, solvents, bleaching or abrasive agents.
- Cleaning materials used should be patient safe and biodegradable.



Wash at 30-74°C. Tumble dry on low heat. Do not iron.
Do not use bleach. Do not use fabric conditioner

DISPOSAL INSTRUCTIONS

- Non-recyclable.
- Dispose of as clinical waste if it has been used by an infectious person.
- Otherwise dispose of through normal waste management.



MANUFACTURER

Hospital Direct (Marketing) Ltd, Units 2, 3 & 4, The Green, Clun, Shropshire SY7 8LG, United Kingdom

EU REPRESENTATIVE: Obelis S.A., Bd Général Wahis 53, B-1030 Brussels, Belgium



www.hospitaldirect.co.uk

HOSPITAL DIRECT (HD) ADDITIONAL SYMBOLS

At HD we understand that not everyone is familiar with all equipment, so to help therapists and healthcare professionals in assessing a product's suitability for a client, we have included two assessment tools, Functional Independence Measure (FIM) and Easy Guide Symbols. These will enable the healthcare professional to assess their client and decide whether or not a product will help the person based on their level of dependence and functional ability.

FUNCTIONAL INDEPENDENCE MEASURE (FIM) ASSESSMENT


Our FIM assessment guide for each product will help you decide the suitability of this product for the person's ability and need. Based on the standard criteria from Level 6 where the person can use the product unaided and unsupervised to Level 1 where all the assistance is provided by the carer and the client can do nothing, this guide is easy, quick and simple to use to check suitability against ability and circumstance.

FIM LEVELS QUICK GUIDE. SUBJECT VS CARER HELP



EASY GUIDE SYMBOLS

We use four symbols to indicate where a product is suitable and safe for the person to use unassisted, with a single carer or with multiple carers. These symbols indicate the minimum recommendation.

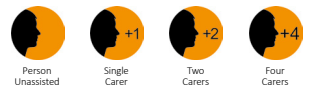

HD has a policy of continuous development and, as such, reserves the right to alter specifications (including measurements, materials and colours) without prior notice. If any serious incident has occurred in relation to the device it should be reported to the manufacturer and the competent authority in which the user and/or patient is established.