



EXPERTS IN PATIENT HANDLING SOLUTIONS



# TRANSFER BOARD SLIDE – WASHABLE

This guide covers all products with the prefix TBS in the product code.

Simple low-cost device that allows for a sitting transfer across boards without risk of shearing induced tissue damage.


## USER INSTRUCTIONS



1

Place board securely on both sitting surfaces, ensuring a safe contact on both.

Place the **Transfer Board Slide – Washable** on top of the board, with the label nearest the leading hip and anti-slip on the board, fabric on the top.

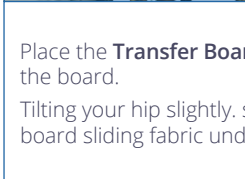

2

Place the **Transfer Board Slide – Washable** lying directly on the board.

Tilting your hip slightly. start to push the top layer of the transfer board sliding fabric under your bottom.



3

Continue to slide the fabric under your bottom.

If you cannot tilt your hip, do not worry, just push as much slide fabric underneath you as possible pressing into the mattress (or cushion if seated) a little to help.



4

Once you have pushed enough fabric (as you judge) underneath your leading buttock, you can begin to start your slide onto the board, using the slide sheet material to help.

Continue sliding so that both buttocks are on the board. You will need to support yourself as shown above in order to make the slide easier.



5

Then, keep sliding across the board onto the other surface, at your chosen, and most comfortable, speed.

You can then let your legs follow and place on the footplate/floor as you would do normally.



6

Once comfortable, remove the **Transfer Board Slide – Washable** by pulling on the slide sheet material, until the product is free. Even if you cannot tilt your hip, the product can still be pulled out easily.

Finally remove the board.



## MATERIAL CONTENT

Fabric: 100% Polyester  
Anti-slip Panel: 100% Polyester with 100% PVC coating



## STORAGE AND HANDLING

- Store in a clean dry environment that is not subject to excessive heat.
- Keep away from sunlight.
- Store flat and avoid placing heavy objects on top of the product.

## WARNINGS AND PRECAUTIONS

- Do not leave or drop on the floor; it is a slip hazard.
- Test for slippiness before use. If over-laundered and doesn't slide please discard or replace.
- Inspect product before use for signs of damage or tears.

## CLEANING INSTRUCTIONS

- Wash as per laundry instructions.
- Can be wiped clean, between use with the same patient, with natural water and detergent or wipes that don't contain alcohol, solvents, bleaching or abrasive agents.
- Cleaning materials used should be patient safe and biodegradable.

## DISPOSAL INSTRUCTIONS

- Non-recyclable.
- Dispose of as clinical waste if it has been used by an infectious person.
- Otherwise dispose of through normal waste management.



Wash at 30-74°C. Tumble dry on low heat. Do not Iron.  
Do not use bleach. Do not use fabric conditioner.



### MANUFACTURER

Hospital Direct (Marketing) Ltd, Units 2, 3 & 4, The Green,  
Clun, Shropshire SY7 8LG, United Kingdom

EU REPRESENTATIVE: Obelis S.A., Bd Général Wahis 53,  
1030 Brussels, Belgium



[www.hospitaldirect.co.uk](http://www.hospitaldirect.co.uk)

## HOSPITAL DIRECT (HD) ADDITIONAL SYMBOLS


At HD we understand that not everyone is familiar with all equipment, so to help therapists and healthcare professionals in assessing a product's suitability for a client, we have included two assessment tools, Functional Independence Measure (FIM) and Easy Guide Symbols. These will enable the healthcare professional to assess their client and decide whether or not a product will help the person based on their level of dependence and functional ability.

### FUNCTIONAL INDEPENDENCE MEASURE (FIM) ASSESSMENT

Our FIM assessment guide for each product will help you decide the suitability of this product for the person's ability and need. Based on the standard criteria from Level 6 where the person can use the product unaided and unsupervised to Level 1 where all the assistance is provided by the carer and the client can do nothing, this guide is easy, quick and simple to use to check suitability against ability and circumstance.

### EASY GUIDE SYMBOLS

We use four symbols to indicate where a product is suitable and safe for the person to use unassisted, with a single carer or with multiple carers. These symbols indicate the minimum recommendation.

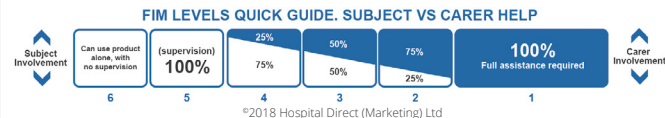

HD has a policy of continuous development and, as such, reserves the right to alter specifications (including measurements, materials and colours) without prior notice. If any serious incident has occurred in relation to the device it should be reported to the manufacturer and the competent authority in which the user and/or patient is established.