



EXPERTS IN PATIENT HANDLING SOLUTIONS



ROTA CUSHION – ESSENTIAL RANGE

This guide covers all products with the prefix RCE in the product code.

The **Rota Cushion – Essential Range** is used for sitting patient rotation and is covered in a cost effective soft breathable fabric that improves comfort and reduces sweating.

USER INSTRUCTIONS



In these instructions we show you how the product is used on a bed, but the same process can be used on a chair or sofa or similar surface that you need to rotate your bottom on.

	<p>1 Place the Rota Cushion – Essential Range on the bed, chair or sofa etc. NOTE: The person must have upper sitting balance, cognitive understanding and some standing ability.</p>
<p>Place your hands next to you if you can, and taking the weight off your seat, allow your bottom to rotate so that your whole body turns.</p>	<p>2</p>
	<p>3 Once turned fully, so you are facing the foot of the bed, you can remove the Rota Cushion – Essential Range.</p>
<p>To remove, rock slightly towards the center of the bed and remove, as shown.</p>	<p>4</p>

MATERIAL CONTENT

Outer Fabric: 100% Polycotton
Filling and Inner Fabric: 100% Polyester
Anti-slip outer: PVC coated polyester, phthalate free
Mini rotating disk: Stainless Steel

STORAGE AND HANDLING

- Store in a clean dry environment that is not subject to excessive heat.
- Keep away from sunlight.
- Store flat and avoid placing heavy objects on top of the product.

WARNINGS AND PRECAUTIONS

- Do not leave the on the floor when not being used.
- Do not use on a wet floor .
- Inspect before use for signs of wear or damage.
- Always test product grip on your chosen surface (bed, car, chair or sofa) before using.

CLEANING INSTRUCTIONS

- Remove the Mini Rotating Disk before washing.
- Cleaning materials used should be patient safe and biodegradable.

DISPOSAL INSTRUCTIONS

- Non-recyclable.
- Dispose of as clinical waste if it has been used by an infectious person.
- Otherwise dispose of through normal waste management.

TOP CUSHION



BASE



Anti-slip Panel



Keep dry



Keep away from sunlight



Caution



Do not use if package is damaged



Storage Temperature Limit



Wash at 30-74°C. Tumble dry on low heat. Do not iron.
Do not use bleach. Do not use fabric conditioner



Patient Information website



Medical Device



MANUFACTURER

Hospital Direct (Marketing) Ltd, Units 2, 3 & 4, The Green, Clun, Shropshire SY7 8LG, United Kingdom

EU REPRESENTATIVE: Obelis S.A., Bd Général Wahis 53, 1030 Brussels, Belgium



www.hospitaldirect.co.uk

HOSPITAL DIRECT (HD) ADDITIONAL SYMBOLS

At HD we understand that not everyone is familiar with all equipment, so to help therapists and healthcare professionals in assessing a product's suitability for a client, we have included two assessment tools, Functional Independence Measure (FIM) and Easy Guide Symbols. These will enable the healthcare professional to assess their client and decide whether or not a product will help the person based on their level of dependence and functional ability.

FUNCTIONAL INDEPENDENCE MEASURE (FIM) ASSESSMENT


Our FIM assessment guide for each product will help you decide the suitability of this product for the person's ability and need. Based on the standard criteria from Level 6 where the person can use the product unaided and unsupervised to Level 1 where all the assistance is provided by the carer and the client can do nothing, this guide is easy, quick and simple to use to check suitability against ability and circumstance.

FIM LEVELS QUICK GUIDE. SUBJECT VS CARER HELP



EASY GUIDE SYMBOLS

We use four symbols to indicate where a product is suitable and safe for the person to use unassisted, with a single carer or with multiple carers. These symbols indicate the minimum recommendation.

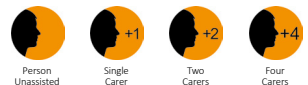

HD has a policy of continuous development and, as such, reserves the right to alter specifications (including measurements, materials and colours) without prior notice. If any serious incident has occurred in relation to the device it should be reported to the manufacturer and the competent authority in which the user and/or patient is established.